



Spice up your workflow

CT motion SPICY




CDadapt™

License

Automatic contrast density reduction (e.g. for cardiac examinations) and simultaneously reduced costs for contrast media due to lower contrast media consumption.
Flexible management of contrast media densities by sequential boli of NaCl and contrast media.



1L GE contrast media Pluspak compatibility plus heat maintaining

Ideal media stock for high patient throughput.
In total 4L (2x 1L contrast media + 2L NaCl).



Enhanced pressure and flow performance

For examinations with high flow rates and high contrast media viscosities.
Increased pressure values up to 20 bar for an increased flow performance.



Extended performance data of RIS/PACS Interface

License

Automatic media label recognition via barcode reading.
Barcode reading of media labels and automatic data transfer to PACS (LOT number, expiry date and GTIN).




Remaining volume indication

Easy-to-understand traffic light indication and clear display of available media volumes plus accurate ml indication.

Transfer of contrast media and NaCl parameters into manufacturer-independent PACS and/or DMS.

By means of DICOM Secondary Capture or DICOM Performed Imaging Agent Administration SR.

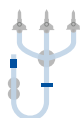

Easier to setup New pump tubing XD 8151

Disposable

Simplified setup procedure due to less insertion steps.
Further improvement of the patient safety due to improved particle filter performance.